


## *400 Series Chemetron Quick-Connect Medical Gas Outlets*

### Elements

- Safety keyed to help prevent cross-connections
- Conforms to NFPA 99
- Rough-in channel box protected by metal plaster guard until final assembly
- Adjustable valve mechanism to compensate for variations in plaster thicknesses up to 3/8"
- Cleaned for oxygen service
- UL Listed

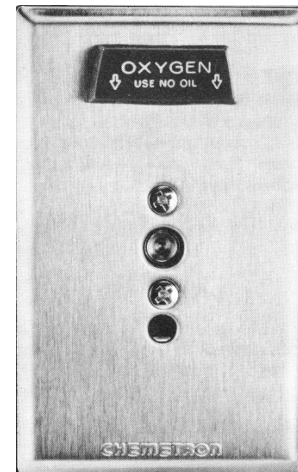

### Description

The Chemetron® 400 Series outlet stations provide convenient access to medical gases in patient care locations. They are fully compatible with Chemetron-style latch-keyed adapters.

Finish plates are stainless steel with a satin finish. Outlet finish plates and rough-in plates are safety keyed for positive cross-connection prevention. Primary valves are separate from the finish plate, allowing easy testing for leaks and cross connection. Primary valves are made of brass, are chrome plated, and contain a self-sealing poppet. Primary valves are field adjustable to compensate for variations in plaster thickness of 1/2" to 3/4".

The type K copper inlet tube has a 3/8" OD connection, and can rotate 360 degrees to

accommodate any field piping arrangement. The primary gas check can be serviced without shutting off the gas supply.

The outlet stations are fully compliant with NFPA 99 1999 and CSA Z-305, and are cleaned for oxygen service in accordance with CGA G-4.1. The outlet is UL listed.

The latch release is color-coded and labeled for the gas service. Release is easily accomplished with one hand using a latch release. The primary valve poppet provides a dust-proof seal. Vacuum outlets may be used with or without slide brackets for vacuum collection bottles of up to 1/2-gallon capacity. (Where slides are used, provide one per vacuum outlet.)



**400 Series Chemetron Quick-Connect Medical Gas Outlets**

**Materials**

<u>Component</u>	<u>Material</u>
Valve Body	Machined Brass Viton Seal
Valve Block	Brass Forging/Machined Brass Type K or L Seamless Copper Tube Viton Seal
Outlet Box	20 Gauge Steel
Cover plate	Stainless Steel #4 Satin Finish Buffed Edge Plastic Color-Coded Release Button

**Ordering Information**

Gas Service	Top Inlet Connection	Quantity	Bottom Inlet Connection	Quantity
Oxygen	64-20-0001		64-26-0001	
Vacuum	64-20-0005		64-26-0005	
Air	64-20-0002		64-26-0002	
Nitrous Oxide	64-20-0006		64-26-0006	
Evacuation Vacuum	64-20-0008		N/A	
Carbon Dioxide	64-20-0004		64-26-0004	
Oxygen/CO2 (CO2 less than 7%)	64-20-0003		64-26-0003	
Slide Bracket	64-59-1212		64-59-1212	





*400 Series Chemetron Quick-Connect Medical Gas Outlets*



**MEDICAL | GAS RESOURCES**  
MEDICAL GAS INNOVATIVE SOLUTIONS

16371 Gothard Street, Suite "F"  
Huntington Beach, CA 92647

**888-743-8298**

[www.medicalgasresources.com](http://www.medicalgasresources.com)  
[info@medicalgasresources.com](mailto:info@medicalgasresources.com)

All specifications are nominal and subject to change without notice.

*Product Specification*